


# RIGHT-FIT BUCKY SYSTEM

## Product Advantages

- Full-time drop protection
- Holder fits into table and wall bucky
- Simple and easy tray rotation (portrait, landscape and top-reference available)
- Optional Removable Grid
- Secure Twin-Lock (trademark) mechanism on holder
- Ergonomic handle on long axis of holder
- Sensors and LED indicators for grid in/out (optional), detector in/out and orientation

## System Components

- Bucky (left or right-hand loading)
- Harmony Table Bucky (left or right head-end)
- Right-fit Mobile DR Panel
- Holder

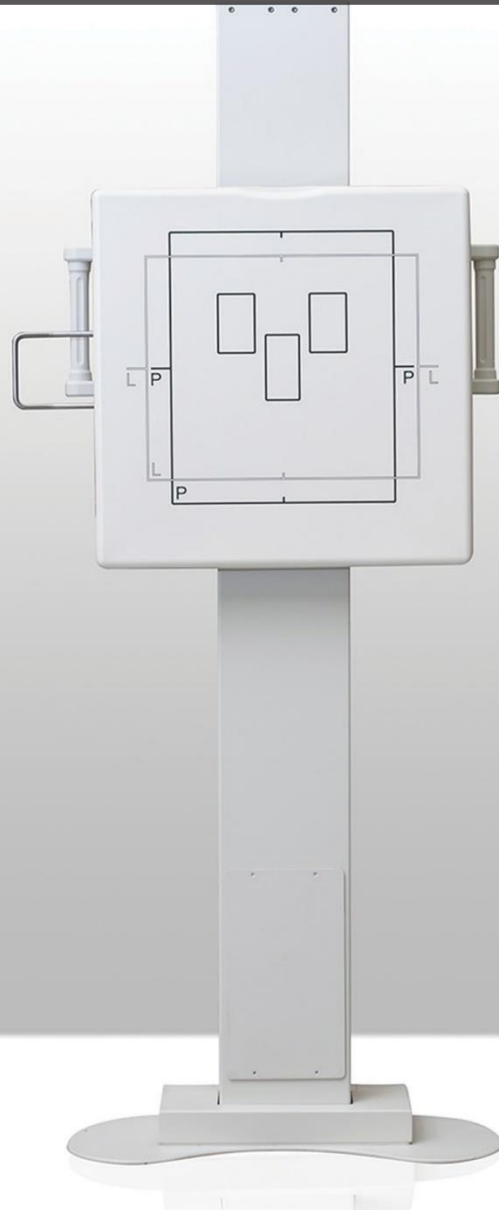

\*\*\*floor column not included\*\*\*

The **Right-Fit Bucky System** is the ideal in-room solution to reduce drops when transitioning the detector between positioning devices or for free exposure use. Your detector will stay safe and secure in the Right-Fit Mobile DR Panel holder, which features an ergonomic handle on the long axis, Twin-Lock (trademark) retaining system and Air-Ride Corners (trademark) for additional drop protection. The Right-Fit Bucky System reduces the TCO (total cost of ownership) of your DR panel by increasing its uptime and useful product lifecycle through damage reduction. It's the only product on the market today that allows full-time use of the Mobile DR Panel Holder, both in and out of the bucky. Take the worry out of a fast-paced radiographic environment and protect your detector with this one-of-a-kind solution for new installations, retrofitting or CR to DR transitioning.

*Additional Notes: Currently available for any 14x17in/35x43cm wireless imaging receptor. Column not included.*

**Reina**Imaging



CASSETTE REPAIR

*Reimagine Radiography*

[www.reinaimaging.com](http://www.reinaimaging.com)





Right-Fit Mobile DRP Holder



Wall Bucky



Table Bucky

## Technical data

	TABLE models	WALL models
Detector rotation	Center rotation	Center / Top-Reference rotation
Loading direction	Left (standard) or Right head end	Left / Right
Flat panel detector tray	Compatible with all ISO 14x17in. wireless image receptors	
Sensors	Detector IN, Portrait/Landscape, Grid IN, Grid type detection (optional)	
Sliding tray grips	Short or Long handle, flat or curved	
Regular power consumption		0.1 A @ 24V
Weight	± 23kg	± 24kg
Size (width x depth x height)	571mm x 584mm x 69mm (Liebel Flarsheim base)	610mm x 610mm x 77mm
Mobile DR Panel Holder	Compatible with all ISO cassette-size wireless image receptors	
Customized models are available upon request		